

GMvCell Help

Overview:

GMvCell simulates single molecule dynamics in a virtual living cell producing sequences of fluorescence light microscopy images built according to the simulated imaging conditions (e.g. illumination method microscope magnification and objective numerical aperture and camera settings). The concept is illustrated below (Fig. 1)

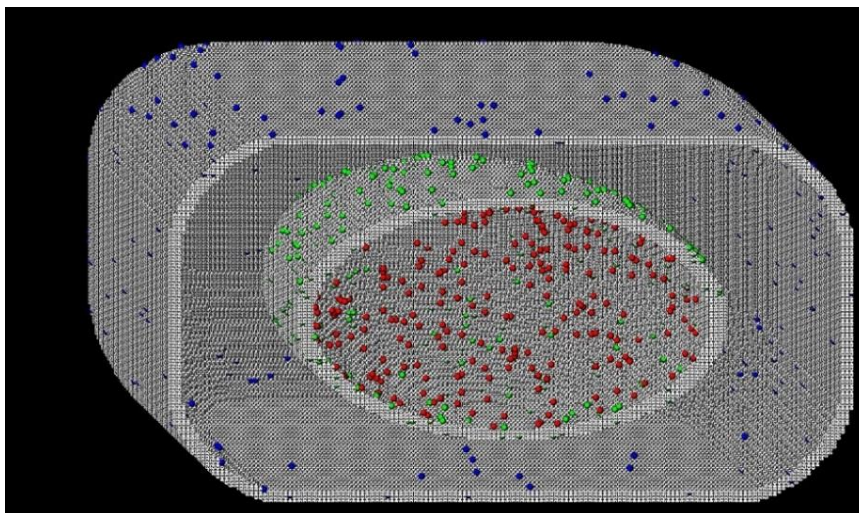


Fig. 1. A voxelated representation of a living cell is generated in computer memory (3D grey outline) and individual molecules (red, green and blue spheres) then move throughout the space in a manner governed by their physical-chemical properties.

The model works with an arbitrary number of molecules belonging to 3 major classes:

- *Cytoplasmic molecules* (“*CytoA*” and “*CytoB*” class): Objects move by a 3-dimensional random walk inside the solution-type voxels. They can bind to and dissociate from objects belonging to “*Memb*” or “*Stat*” classes. If binding occurs the pair move together with mobility determined by the slowest object (usually in membrane-type voxels).
- *Membrane localised molecules* (“*Memb*” class (*A*, *B*, and *C*)): Objects are localised in the membrane-type voxels and can move by a 3-dimensional random walk (in correct type voxels) or can be immobile. They can bind/unbind *Cyto* and other *Memb*-type class objects.
- *Static molecules* (“*StatA*” class): Objects usually located inside the cell as groups of immobile molecules forming fibres or filaments (e.g. cytoskeletal components). They can bind/unbind *Cyto* class objects.

Each object may have up to 4 green and 4 red fluorophores. At each time point (virtual video frame) all object positions are computed according to a Brownian random walk, calculated knowing the properties of each object. Object collisions are recorded and binding/unbinding events are computed

based on stochastic first order kinetics. The new object positions are then used to compute the number of photons emitted by each fluorophore, depending on the level of illumination at that particular point in the cell. The emitted photons are collected by a virtual objective that produces a diffraction-limited point spread-function at the virtual camera. Simple geometric optics are computed so that a virtual image of the cell is formed by the emission from all active (non-bleached) fluorophores illuminated under the chosen imaging conditions (including counts due to additional noise). The sequences of images, generated during the model run, can be saved for future analysis as a set of BMP files or as a single GMV movie – a proprietary file format used in the partner software for single molecule detection, tracking and export - “GMimPro” www.mashanov.uk or <https://github.com/GMashanov>

Detailed information about this model is published in *Mashanov and Molloy (2024)* <https://www.nature.com/articles/s41598-024-70925-2.pdf>.

The description of modelling single molecule fluorescence and binding kinetics is published in Mashanov (2014) [Single molecule dynamics in a virtual cell: a three-dimensional model that produces simulated fluorescence video-imaging data - PMC \(nih.gov\)](https://pubmed.ncbi.nlm.nih.gov/25111111/)

The program flow-chart is shown in Fig. 2 (below).

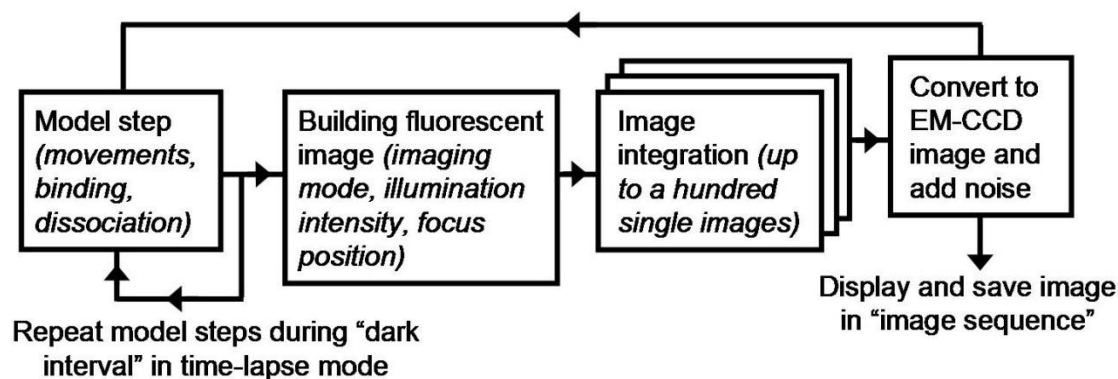


Fig. 2. Block diagram showing the computation cycle used to generate simulated fluorescence microscopy movies of single molecules moving within the cell volume (see Fig. 1).


In the user guide below the following colour codes are used:

- **BLUE UPPERCASE** indicates top-level menu tabs
- **Red** is used for sub-level menu tabs.
- ***Black bold italic*** font is used for parameter entry headings
- **Magenta** is used for short-cut key options


Running the program:

The cell structure and molecules must first be initialised via the **SETTINGS** tab on the main screen or accessed using the **F2** function key.

The **SETTINGS** dialog window (also accessed via **F2**) contains several tabbed pages that are used to set the properties of each class of objects the virtual cell parameters and the simulated imaging conditions.

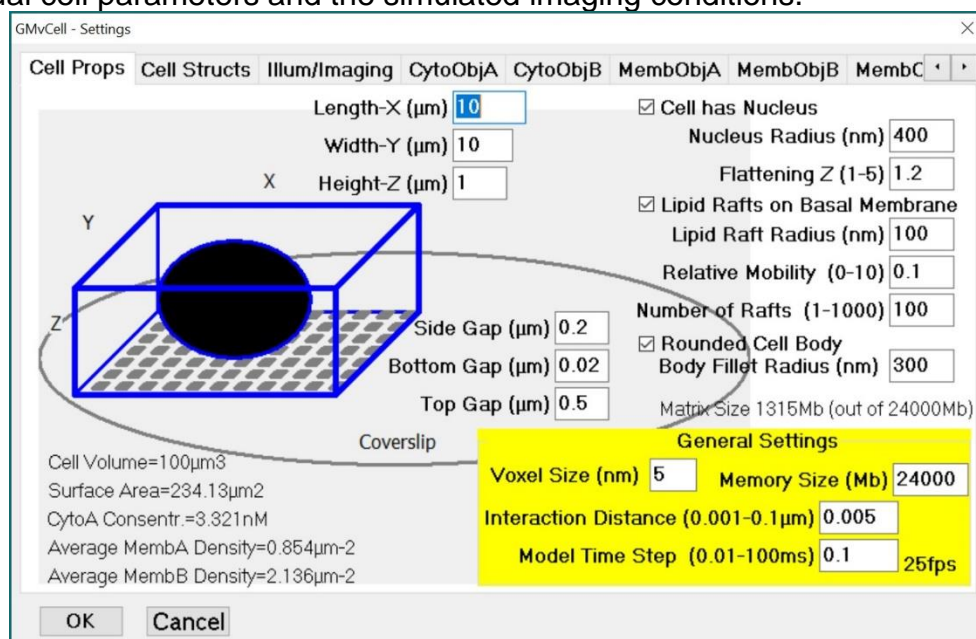

Cell Props: The virtual cell can be a rectangular cuboid of a set cell size of **Length-X**, **Width-Y**, and **Height-Z** dimensions (in μm). The lower or basal surface is composed of plasma membrane, is closest to the microscope objective and is set to be at the focal plane (defined as $\zeta=0$). Cell is placed in voxelated 3D matrix with outside gaps defined by corresponding **Side Gap**, **Bottom Gap** and **Top Gap** values. The coordinate system is mapped to computer memory as described in the *Main Paper* and is shown graphically below:

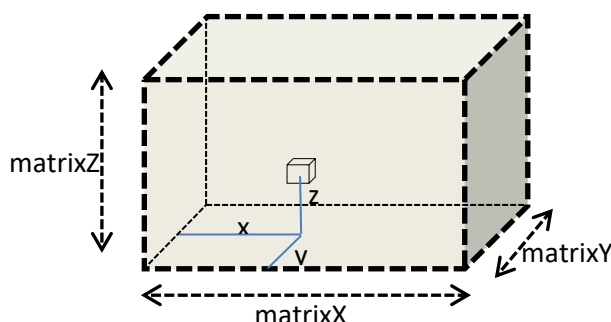

The cell may optionally include an ellipsoidal nucleus **Cell has Nucleus** with a **Nucleus Radius** (in nm) in x-y plane and an extended or flattened radius in the z plane, **Flattening Z** (where values <1 give an oblate spheroid and >1

prolate; range is 0.05 to 2). The user must ensure that the nucleus size (in all dimensions) does not exceed the cell size.


The cell membrane may contain circular **Lipid Rafts on Basal Membrane** with a **Raft Radius** (in nm) Membrane associated protein molecules moving in these voxels either slow-down (<1) or speed-up (>1) according to the **Relative Mobility** value. User must set **Number of Rafts**, which will be initiated in MATRIX submenu.

The cell body may optionally have rounded edges by selecting the check box option **Rounded Cell Body** with a **Body Fillet Radius** entered in nanometres.

General Settings found in the lower right, yellow highlight region are used to set the **voxel size (nm)** and allocated **memory size** in megabytes (you need to close the "Settings" window by clicking "OK" to set this value). **Interaction Distance** (i.e. average size of the interacting objects (molecules), given in μm (e.g. 3-8 nm for protein molecules) for general calculations – bigger the sizes - bigger the probability of binding (also depends on mobility values).

Model Time Step is the time interval used for each calculation cycle, Δt . Many time-steps can be summed to build one fluorescent image/frame (see Illum/Imaging Tab). Note that if the model simulates fast moving molecules (e.g., diffusion rate $2\text{-}20 \mu\text{m}^2 \text{s}^{-1}$) and frame rate is slow (<50 fps) it is advised to sum a few short time-steps (e.g., sum of 100 time-steps 0.5 ms each) to create one fluorescent image (see Illumination/Imaging page [N of Steps summed in One Image](#)).

In the left corner of **Cell Props** tab you can see some statistics: Approximate cell volume and surface area (all 6 sides of the cell), concentration of CytoA objects, and density of membrane (A/B) objects.

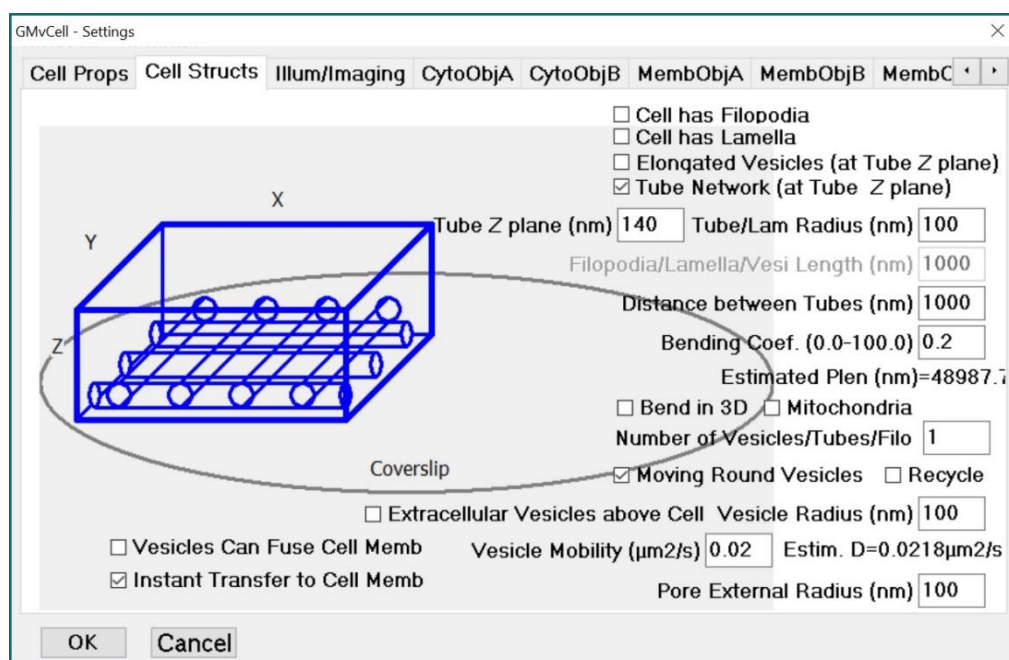

Cell Structs: The virtual cell can have filopodia or a body-wide lamella protruding from the cell on the right side at bottom level (simulating a lamellipodium). Filopodium diameter or lamella height are set in the **Tube/Lamella Rad. and Length** in **Filopodia/Lamella/Tube Length**. The number of filopodia is set in **Number of moving Vesicles** edit window. A number of elongated vesicles (enclosed tubes) separated by **Distance between Tubes** can be constructed inside the cell oriented from left to right, tilted upwards according to **Tilting/Bending Coefficient**. A randomly oriented interconnected tube network can be constructed in 2D (at $Z = 2 \cdot \text{tubeRad} + \text{tubeGap}$ height) or in 3D. Tubes will randomly bend according to the **Bending Coef.** value. In addition, **Moving Round Vesicles** are placed either inside the cell body or above the body in case of **Extracellular Vesicles** (make sure that **Top Gap** value is bigger **Vesicle Radius * 2**). **Vesicle Mobility** can be set manually, the estimated default mobility, D is shown for a sphere of given radius, moving in water at 25°C. **Pore External Diameter** is used for Exocytosis simulation when vesicles may fuse to the plasma membrane (either from the inside or outside), except for the situation where **Instant Transfer to Cell Membrane** is selected, in which the vesicular membrane molecules are instantly moved to the nearest plasma membrane voxels and molecules within the vesicle lumen move to the nearest extracellular/cytoplasm voxels (depending from which side the vesicle is fusing). In case of Instant Transfer, empty vesicles are removed from the matrix and replaced by cytoplasmic or extra-cellular voxels. **NOTE: Moving Round Vesicles** mode cannot be run together with **Elongated Vesicles** because they share the same properties.

GMvCell - Settings

Cell Props Cell Structs Illum/Imaging **CytoObjA** CytoObjB MembObjA MembObjB MembObjC

Number of CytoA Objects CytoA Concentr.=110.703nM, 66.667 molecules/ μm^3

N of Green Fluorophores per CytoA N of Red Fluorophores per CytoA

Green Fluorescent Fraction (%) Red Fluorescent Fraction (%)

Diffusion Rate ($\mu\text{m}^2/\text{s}$) Aver. Jump (RMSD) per TmStep=122.47nm

Initial Distribution of CytoA

☒ CytoA polymerase into static filaments

CytoA Polymerisation Rate (s-1) CytoA Depolymerisation Rate (s-1)

OK Cancel

CytoObjA: This menu tab is used to set the properties of objects moving within the cytoplasm: **Number** of CytoA objects (up to 50000 objects),

number of **Green Fluorophores** per CytoA object and **Fraction** of CytoA objects with active fluorescent tag. Note if an object has more than one green fluorophore per molecule and fraction is set to <100%, some CytoA objects may have a different number of active fluorescent tags at the beginning of the model run. The same rule applies to the **Red Fluorophores** on CytoA objects (and to the **MembA&B&C** and **StatA** objects, see below). The user cannot alter the number of objects or fluorescent tags number in the middle of a model run (the model will restart) but can change the **Diffusion Rate**. CytoA objects are placed at random positions in voxel matching type of Initial Distribution, which can be cytoplasm, nucleoplasm, tubes or vesicles lumen, and in extracellular space. **Note:** The user must ensure that structures are reset before running the model and click **Reset** on the main screen. CytoA objects can polymerize into static filaments using **Polymerize** **Depolymerize** rate constants.

CytoB: This object class is similar to CytoA (see above for detailed description) (up to 50000 objects), but have fewer Initial distribution options (cytoplasm, outside cell, and in cell centre). CytoB can bind/dissociate from CytoA molecules according to Bind / Dissociate rate constants. Fluorescent CytoB can be added to simulate background noise that adds in addition to sources of camera noise (see below).

GMvCell - Settings

Illum/Imaging CytoObjA CytoObjB **MembObjA** MembObjB MembObjC StaticObjA

Number of MembA objects: 2000 Average MembA Density=34.069 μ m⁻²

N of Green Fluorophores per MembA: 1 N of Red Fluorophores per MembA: 0

Fluorescent Fraction (%): 100 Fluorescent Fraction (%): 100

MembA Diffusion Rate (μ m²/s): 0.1 Aver. Jump (RMSD) per TmStep=54.772nm

Initial Distribution of Memb A: Random - Cell Membrane

☐ MembA can bind CytoA ☒ MembA is an Organelle Pump

CytoA-MembA Binding Rate (s⁻¹): 1000000 CytoA-MembA Dissociation Rate (s⁻¹): 0.2

☒ MembA can polymerise

MembA-MembA Binding Rate (s⁻¹): 100000 Max MembA Units in Polymer: 8

MembA-MembA Dissociation Rate (s⁻¹): 1

OK Cancel

MembObjA / B / C: These three pages contain settings analogous to the CytoA page (see above) (up to 50000 objects), but its **Initial Distribution** ComboBox has selection corresponding to membrane structures (cell membrane, nuclear membrane, tube, and vesicle membranes). There are also tick-boxes which control binding of Cyto and Memb molecules to each other. MembA molecules also can polymerize (up to 16 monomers). The binding and dissociation events are determined by **Binding Rate** and

Dissociation Rate coefficients. For the binding events, the model checks the Average Interaction Time (AIT) for every eligible pair of molecules. AIT depends on the distance between molecules, their mobility, and **Interaction Distance** (set on CellProps page). The AIT is multiplied by Binding Rate and used as an input for RNG determining the outcome of binding test. The **Dissociation Rate** is used as a parameter in RNG testing possible dissociation event for every bound pair of molecules every time step.

GMvCell - Settings

Illum/Imaging CytoObjA CytoObjB MembObjA MembObjB MembObjC StaticObjA

Number of StatA objects 12000

Green Fluorophores per StatA 2 Red Fluorophores per StatA 0

Fluorescent Fraction (%) 10 Fluorescent Fraction (%) 100

Distribution of StatA Qbjects Straight Fibers in XY (at z=100nm)

Distance between Objects in a Fiber (nm) 8

Monomers Per Fiber 200 Fiber Extension (nm) 50

☒ Static A can bind Cyto A

CytoA-StatA Binding Rate (s-1) 500000 CytoA-StatA Dissociation Rate (s-1) 0.1

☒ CytoA Can Move to the right along StatA

Motor Step Size (nm) 8 Motor Stepping Rate (s-1) 100

OK Cancel

StaticObjA: Static objects remain stationary and cannot change position inside cell during the model run (up to 100000 objects). The “StatA” objects can form “Straight Fibres”, and “Curved Fibres” (laying in horizontal plane at $\zeta=100\text{nm}$), “Fibres Curved in 3D”, and specific case of filament bundles in filopodia (6 monomers place at the center of filopodia core (see **CellStructs** tab). The length of the fibers depend on the **Distance between Objects in fiber** and **Monomers per Filament** value. The length of the filament will be limited if it grows into an obstacle (membrane structure). CytoA molecules bound to StatA filaments can move to the left along the fibres in case of actin bundles and random directions in case of fibers in cell body. Both tick-boxes **CytoA can Bind StatA** and **CytoA can Move to the right along StatA** should be checked. The velocity of CytoA movement determined by **Motor Step Size** and **Motor Step Rate**.

Illumination/Imaging: This page contains a group of parameters. First group controls the emission of single fluorophores. **Diffraction Limited Spot Size** (Full Width at Half Maximum, FWHM) of a single molecule image placed at the focal plane of an objective lens can be set close to the theoretical limit for given wavelength ($\sim 1/2$ of the wavelength of emitted light). There is also **Emission Rates** and **Photobleaching Rates** for both Green and Red fluorophores. NOTE - these rates set for the fluorophores placed in the centre of the image at coordinates = $(x/2, y/2, \zeta=0)$; where illumination is brightest. Fluorophores placed in other parts of the cell could receive much weaker illumination (depends on the illumination pattern) and as consequence emit proportionally smaller number of photons. The bleaching rate for Green and Red fluorophores is also set for the molecules in the centre of the image $(x/2, y/2, \zeta=0)$. It will be proportionally smaller for the weakly illuminated molecules. The image of single fluorophore (2-D Gaussian spot) placed above or below current focal plane will be increased (blurred), while the number of emitted photons remains the same. Ideally the relation between fluorophore FWHM and distance from the focal plane is described by simple hyperbolic function $FWHM(z) = K_{obj} * (\zeta - \zeta_{focal})^2 + FWHM_{focal}$, where ζ is the z-position of a molecule (above coverslip level), ζ_{focal} is a current position of focal plane of the objective lens (usually ζ_{focus} is defined as 0), $FWHM_{focal}$ is the size of a fluorophore image at the focal plane (see above). The measured relation between FWHM of single fluorescent molecule (Cy3B) and movements along the z-axis of the objective lens 100X1.45NA resulted in **Kobj** $\approx 0.6 \mu m^{-1}$.


Illumination: There are 3 modes of illumination used in this model: Epi-illumination, Total Internal Reflection Fluorescence Microscopy (TIRFM), and confocal imaging. It is assumed that in two first modes laser beam has a 2D Gaussian shape profile determined by its FWHM, The user can set **Laser Beam FWHM** to control the pattern of illumination. The FWHM can be increased (e.g., \gg Cell Length and Width) to make the nearly-flat illumination

pattern in x-y plane. The illuminated area can be shaped by a circular **Iris Diaphragm**. In “Epi-illum” and “Confocal” modes the whole cell volume will be illuminated within the limits of the size of **Iris Diaphragm**. The cell profile (along cell length) is shown next to these Edit Controls. The TIR angle parameter can be changed between 63 and 75 degrees to control the depth of TIR illumination. NOTE – It also depends on the wavelength of the illumination (488nm for Green fluorescence and 561 nm for Red fluorescence simulation).


Camera Settings: Few model steps can be summed to create one camera image. It gives realistic images of the fluorescent objects moving considerable distances during one image acquisition. For example, if we set model time step (see above) to 5 ms and set N of **Model Steps summed in One Image** to 10 we will get camera rate of 20 fps which is achievable by many sensitive cameras. At this rate slow moving molecules will look like sharp spots of light, while fast randomly moving molecules look like dim clouds. If the **Number of Dark Time Steps** is set above zero, the model will run in a time-lapse mode when one illuminated image is followed by specified number of dark time-steps. For example, if we use above settings (summing 10 images) and set **N of Dark Time-Steps** to 90, the model will display (and record) one fluorescent image (50ms exposure time) every 0.5 s giving 2fps (The photobleaching rate will be reduced accordingly). The model will add few “noise” counts to every pixel on the image which simulate camera and photon noise. The noise counts have Gaussian random distribution which sigma value is set as a **Camera Noise**. EMCCD cameras multiply every registered photon by an average number determined by **Camera Gain**. This multiplier has some degree of variation which, in this model, was set as a square root of pixel value. The EMCCD mode can be switched ON/OFF using **EMCCD camera** check-box. If two colours are used in the model run the results can be displayed / recorded in few modes – **Two Camera** mode used to display Green and Red images side by side (NOTE - make sure that both lasers are switched ON in the main menu). **Alternating Illumination** mode is used to illuminate one frame with one colour and next one with another (Green-Red-Green-Red and so on). The results can be displayed and recorded as one image (.gmv) or as double (.gmG and .gmR) if **Two Camera** box is checked.

Tip: In order to visualise the above structures during the model run the user should run the model for a few tens of frames to allow “Cyto” molecules to equilibrate across the cell volume then by pressing the “C” key “Cyto” objects are marked as white dots. It is recommended to use a short time steps (e.g. 1 ms) to model “CytoA” dynamics within the tubular structures and around the nucleus because molecular movements are limited in the narrow spaces. It is also advised to limit the cell height (z) so that the tubular network can be illuminated by the evanescent wave. By pressing “C” again the white dots are hidden.

The Main Window menus.

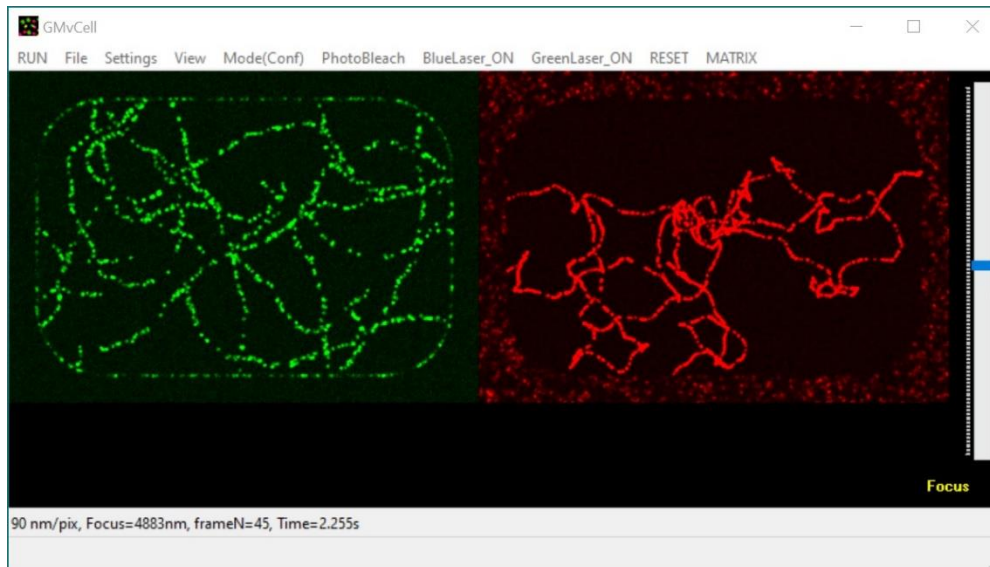

Run submenu: User can start and stop model run but pressing “**RUN/STOP**” menu item or by using “**F9**” / “**Esc**” keys.

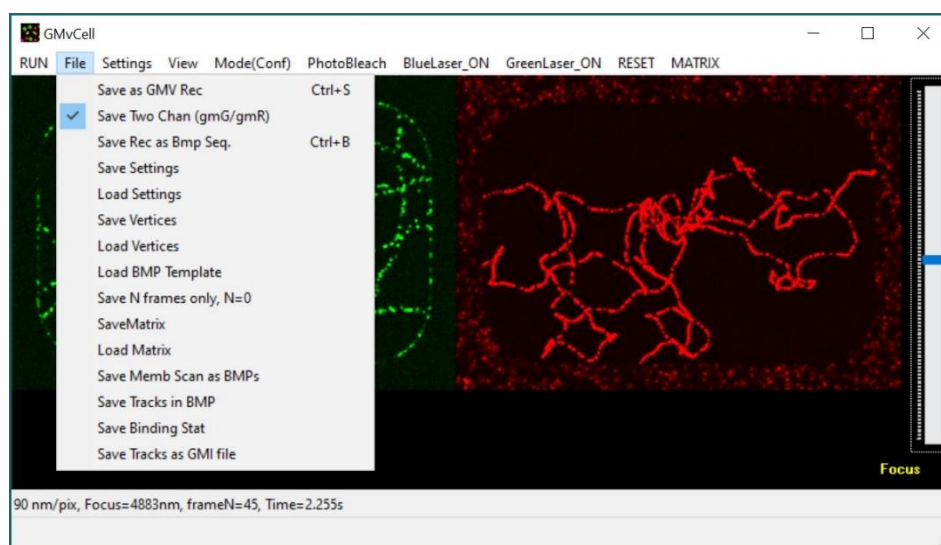

File submenu:

at any moment before or during model run user can start recording “sequence of fluorescent images” in the specific GMV file format used in Single Molecule Detection and Tracking software GMimPro (see separate GMimPro help files and recommended papers to learn this software). GMimPro can be used to display, analyse, and save as sequences of BMP or binary files original GMV files generated by the model. GMimPro webpage (www.mashanov.uk) also contains ImageJ plug-ins for importing GMV files into ImageJ. In case of two colour **Two Camera** recording, image sequences can be saved in two separate files **Save Two Chan (gmG/gmR)**. The results also can be saved as 24-bit BMP files using **Save Rec. as BMP Seq.** Please specify the file name for the first frame. WARNING – file name should not have number or “.bmp”

string at the end of the filename because the index specifying image number will be added to the filename automatically. It is strongly advised to create a separate folder for each BMP image sequence created at the model run. In both record modes the recording is terminated by pressing **ESC** key. WARNING – make sure that the fluorescence (**'F'** key) and at least one of the lasers is switched ON before you start recording.

Most initial conditions in the **Settings** menu can be stored for future model runs using **Save Settings** option. The current settings will be stored in a text file with default extension “gms”. This is text type file containing text strings encoding model/object properties in following format: *property,value;* (e.g., *cellSizeXnm,5850; raftMobilityRatio,0.02; cytoAdistribution,NUCLEOPLASM;*), where property is separated from value by ‘,’ and terminated by ‘;’. These




settings can be loaded using **Load Settings** line in the File menu. NOTE – if the numbers of objects of any class or cell sizes loaded from settings file do not match the current settings, the model will start a new run using initial distributions of membrane and static objects stored in loaded file. If a set of vertices was created in MATRIX menu (**Set Membrane Vertices**), it can be saved as a file (**Save Vertices**) for future download (**Load BMP template**). A template **8-bit** BMP file, containing 2D projection of 3D membrane structure can be stored as a template of membrane shell, which would have the Z-height of **CellHeight**. The bmp x and y sizes will be stretched/shrink according to **CellLength** and **CellWidth**.


The number of recorded frames can be limited to N, if it is set above 0 (**Save N frames only**). Constructed matrix can be saved to disk and loaded back to model (**Save Matrix / Load Matrix**) to save time required to construct large size structures. The operator has to check that cell and surrounding sizes are the same when matrix is saved / loaded into program. Non-destructive compression is used to reduce the size of the matrix on disk. **Save Tracks as GMI file** - xyz trajectories containing exact positions of MembA objects can be saved in GMI format to be open in “Motility” software, (www.mashanov.uk) to be analysed or exported as txt/CSV data for future analysis.

View submenu

The user can change camera pixel sizes used in the model (default value 100

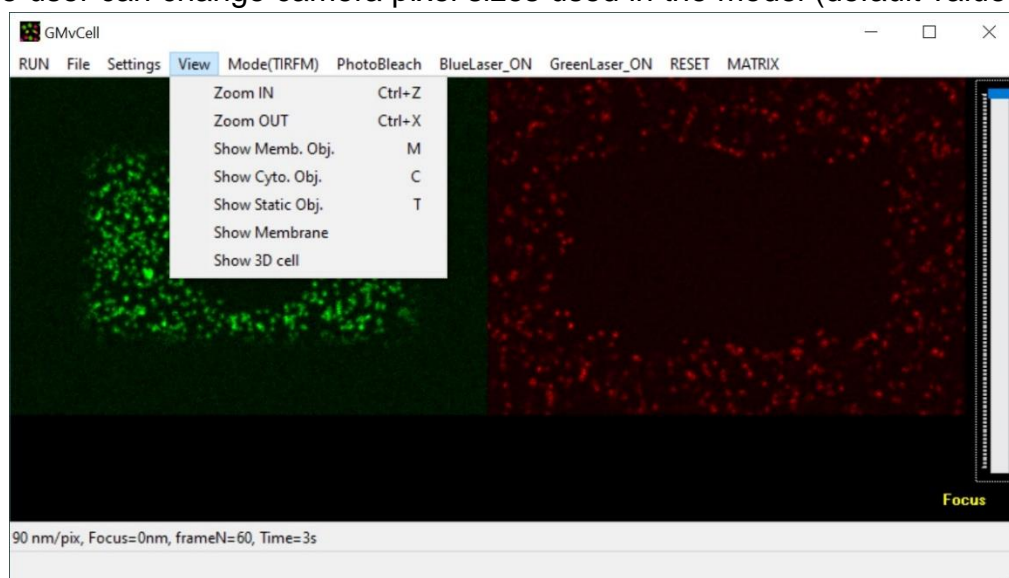

nm/pixel) using “**Ctrl+Z**” and “**Ctrl+X**” key combinations (current scale is shown in status bar at the bottom of main window). NOTE – These changes will affect the sizes and brightness of the images of fluorescent objects projected on camera CCD. Normally, objective lens x100 would increase image size 100 times to the size of CCD pixels (e.g., Andor iXon897BV has CCD pixel sizes 16 μm , so an objective lens x100 would give camera scale 160 nm/pixel). The current XY positions of “Cyto” objects can be shown on CCD image using ‘**C**’ key (white dots), “memb” objects using ‘**M**’ key (red/yellow dots), “stat” objects using ‘**T**’ key (blue dots). You can see membrane structures at current Z plane, if you select **Show Membrane** in a View menu – use “**grey**” “**+**” / “**-**” keys to move Z-plane up and down. Make sure that model is not “running” otherwise membrane image will be overwritten by fluorescent image. You can see 3D cubical representation of membrane structures and molecular objects in a separate “Cell3D” window.

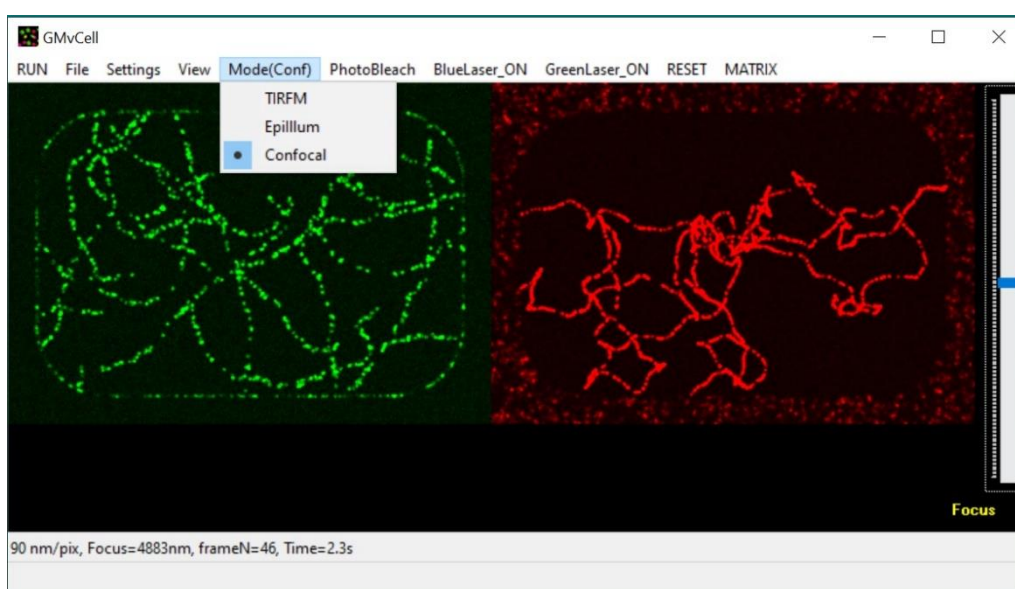


Imaging Mode submenu

TIRFM (Total Internal Reflection Fluorescence Microscopy) – cell is illuminated according to the rules of optics for Evanesce Wave propagation at the interface between mediums with high and low refractive indices (See Axelrod et. al., 1992). The profile of illumination is determined by the colour of illumination, TIR angle, and laser beam profile settings (see Settings->Illum/Imaging menu).

Epillum. – The laser beam penetrates the whole cell. Its XY intensity profile depends on the **Laser Beam FWHM** and **Iris Diaphragm** settings (see Settings Menu).

Confocal mode – It is a simulation of scanning method of illumination / imaging where intensity of illumination is constant at any XY point in the current focal plane (laser beam focused in the focal plane of the objective lens). The intensity of illumination decreases dramatically above and below focal plane (outside the focused laser beam) due to a sharp increase in the size of laser beam cone. It is also assumed that a “virtual pinhole” blocks most of the emission from the fluorescent objects placed above or below current focal plane.

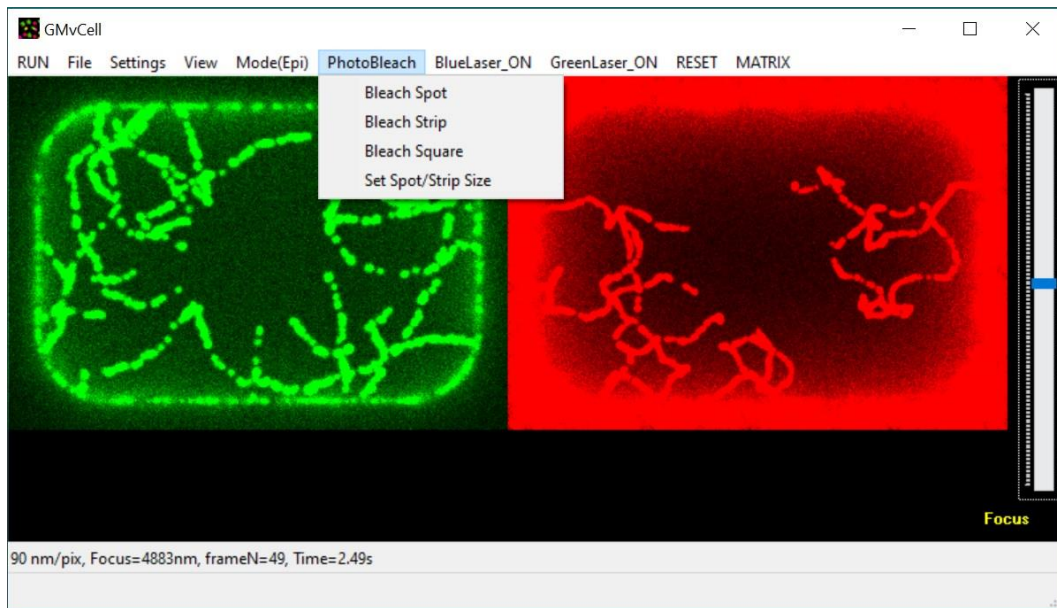


Photobleach: The user can bleach all the fluorescent molecules within chosen spot (circle), square, or strip area. NOTE - If bleaching is done in TIRFM mode “Cyto” and “Stat” molecules will not be affected by the bleach (only molecules at basal cell membrane (“Memb” objects)), but in the “Epillum” or “Confocal” modes molecules of all types will be bleached. Make sure that the one or both lasers (Blue, Green) are ON at the moment of bleaching to bleach one or both colours. The size/diameter of bleaching area can be set in [Set Spot/Strip Size](#) Edit control. The Round Spot or Square is placed in the middle of the cell, while bleached Strip is placed at the right edge of the cell.

BlueLaser and **GreenLaser** menu items allow to switch Blue or Green illumination ON and OFF. These controls will not work if [Alternating Illumination](#) check-box is checked in the [Settings](#) menu. NOTE – make sure that Illumination is ON (key ‘**F**’) to bleach molecules.

RESET - click this menu item to reset the model – all the objects will be placed according to the initial distributions and all available fluorophores (see Settings) will be reactivated (made fluorescent again if they were bleached in the previous run).

MATRIX menu: From here you can create the cell structures in a matrix space ([Make Geom. Shape](#)) before user resets (**RESET**) all the molecule objects in 3D space and model run (‘**F9**’). This can take substantial time, if the size of the matrix is bigger 1 Giga voxel (see general Settings on CellProps Tab). Cell structures can be made of 2D .bmp template loaded in File menu. You need to load such file, check [Make shell from BMP](#) line and Make Geom. Shape out of it – it will extrude 3D shape out of 2D shape and hollow it to create membrane structures.

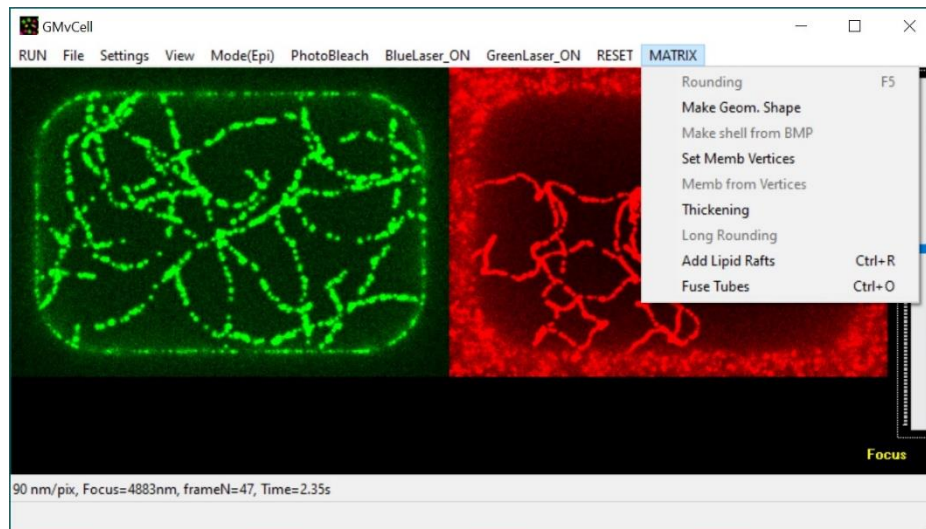

Membrane vertices can be set out of membrane voxels (connected to nearest membrane voxels). The position of these vertices can be rounded ([Rounding](#)) to smooth sharp membrane corners/edges in the matrix. Then membrane can be reconstructed in matrix space ([Membrane from Vertices](#)) with new coordinates. It is important to thicken membrane when you use rounded structures (e.g., tube network, rounded cell, filopodia and so on) to avoid membrane voxels connected to other voxels only at edges to prevent solution molecules leak though the membrane and ensure free movements of membrane molecules.

[Add Lipid Rafts](#) – will add N number of round membrane rafts (see CellProps Tab) at random position to basal cell membrane by changing membrane voxel properties. Membrane molecules, moving in these voxels, will slow down according to the ratio set in CellProp Tab.

List of Key Combinations

/ Shortcuts Start/Stop/Reset:

1. F2 – Launch Settings dialog Window.
2. F4 – Reset time (t=0s)
3. F9 – Run
4. Esc – stop run / stop recording
5. 'Space' make one model step. Note: more than one step may be required to form a new image on the screen.

Action Controls:

6. Ctrl+'S' - Save GMV record
7. Ctrl+'B' – Save record as Sequence of BMP files
8. Ctrl+'Z' – Zoom In – Decrease size of CCD pixel
9. Ctrl+'X' – Zoom Out – Increase size of CCD pixel
10. Ctrl+'N' – Reset the model – go to initial state and distribution of present objects.
11. 'F' – Switch ON/OFF illumination. Note: make sure that one or both lasers are ON to be able to illuminate “the cell”.

Display Modifiers:

12. 'C' – Show/Hide white dots marking x-y positions of CytoA objects
13. 'M' – Show/Hide red and yellow dots marking x-y positions of MembA

and MembB objects

14. 'T' - Show/Hide blue dots marking x-y positions of StatA objects